


How to install the suprafab-analytics mount and unmount buttons

1. Go to the measurement computer you want to install the buttons.
2. Log in into windows using the win admin account.
3. Use the zellbio-users mount button and log in to the zellbio-users network drive.
4. Copy the folder "bcpfs-mount" from the zellbio-users network drive to C:\\
Note: Do not delete the "bcpfs-mount" folder from the zellbio-users drive.
5. Copy the "bcpfs-mount" folder to C:\\
6. Open the Command Prompt (cmd) as admin by pressing the Windows key, type "cmd" and select in the list "Run as administrator" and click "yes".
7. Go in the cmd to the "bcpfs-mount" folder by typing "**cd c:\bcpfs-mount**" and confirm by pressing enter
8. There are two options depending on the win system:
 - a. for Win 10 and higher type in the cmd "**handle-4-1.exe bcp.storage01**" and confirm by enter and in the appearing window select "agree"
 - b. for Win Vista/XP and lower type in the cmd "**handle-3-51.exe bcp.storage01**" and confirm by enter and in the appearing window select "agree"
9. Close the cmd and in case the win user account is different from the win admin account, log off as admin and log in as user.



10. Right click on an empty space on the desktop and select "new" and then "shortcut".
11. Use the Browse button to navigate
 - a. for Win 10 and higher to "C:\bcpfs-mount\mount-suprafab-analytics"
 - b. for Win Vista/XP and lower to "C:\bcpfs-mount\mount-suprafab-analytics-winxp"
12. Click Next and then Finish.
13. Use the Browse button to navigate
 - a. for Win 10 and higher to "C:\bcpfs-mount\unmount-suprafab-analytics"
 - b. for Win Vista/XP and lower to "C:\bcpfs-mount\unmount-suprafab-analytics-winxp"
14. Click Next and then Finish.

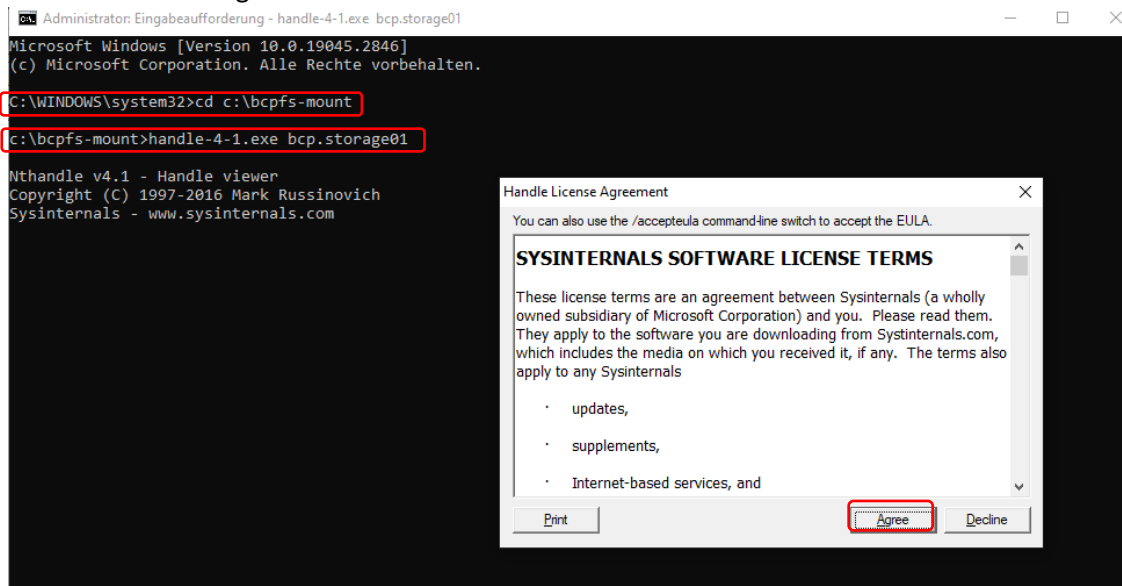