


How to install the suprafab-analytics mount and unmount buttons

1. On your personal computer, download the folder bcpfs-mount from <https://box.fu-berlin.de/s/XGaKn3exnJgMmPZ> and copy it to an empty and formatted USB device
2. Go to the measurement computer you want to install the buttons.
3. Log in into windows using the win admin account.
4. Connect the USB device with only containing the bcpfs-mount folder and copy the folder "bcpfs-mount" to C:\
5. Open the Command Prompt (cmd) as admin by pressing the Windows key, type "cmd" and select in the list "Run as administrator" and click "yes".
6. Go in the cmd to the "bcpfs-mount" folder by typing "**cd c:\bcpfs-mount**" and confirm by pressing enter
7. There are two options depending on the win system:
 - a. for Win 10 and higher type in the cmd "**handle-4-1.exe bcp.storage01**" and confirm by enter and in the appearing window select "agree"
 - b. for Win Vista/XP and lower type in the cmd "**handle-3-51.exe bcp.storage01**" and confirm by enter and in the appearing window select "agree"
8. Close the cmd and in case the win user account is different from the win admin account, log off as admin and log in as user.



9. Right click on an empty space on the desktop and select "new" and then "shortcut".
10. Use the Browse button to navigate
 - a. for Win 10 and higher to "**C:\bcpfs-mount\mount-suprafab-analytics**"
 - b. for Win Vista/XP and lower to "**C:\bcpfs-mount\mount-suprafab-analytics-winxp**"
11. Click Next and then Finish.
12. Use the Browse button to navigate
 - a. for Win 10 and higher to "**C:\bcpfs-mount\unmount-suprafab-analytics**"
 - b. for Win Vista/XP and lower to "**C:\bcpfs-mount\unmount-suprafab-analytics-winxp**"
13. Click Next and then Finish.

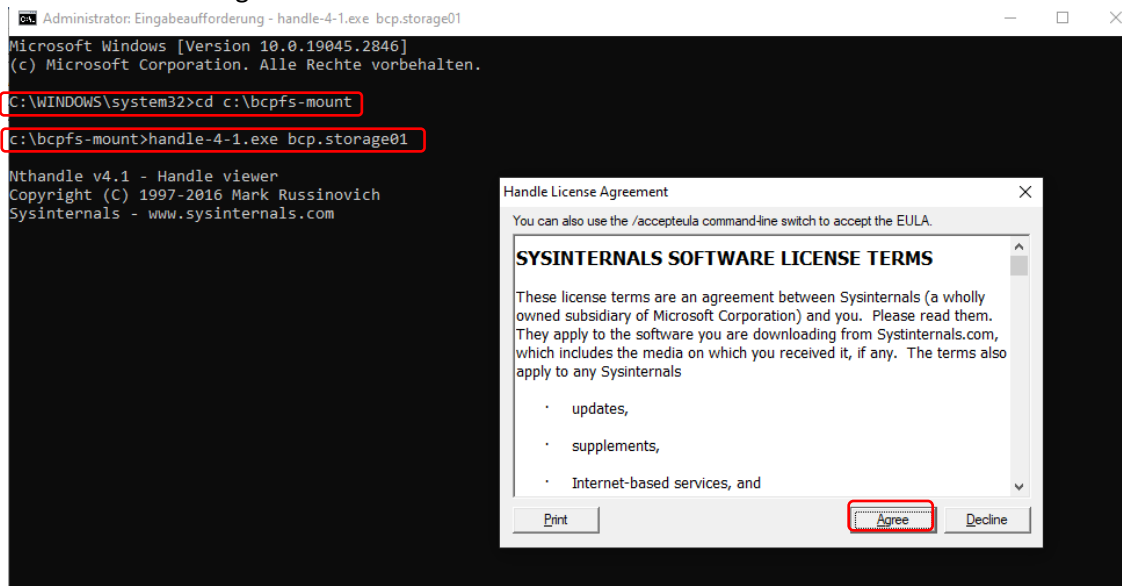